



AMH Rapid Test Cassette (Whole Blood/Serum/Plasma)

Package Insert

REF FAMH-402 English

A rapid test for the semi-quantitative detection of the anti-Müllerian Hormone (AMH) in whole blood, serum or plasma.

For professional *in vitro* diagnostic use only.

[INTENDED USE]

The AMH Rapid Test Cassette (Whole Blood/Serum/Plasma) is a rapid chromatographic immunoassay for the semi-quantitative detection of the anti-Müllerian Hormone in whole blood, serum or plasma to aid in the diagnosis of ovarian functions.

[SUMMARY]

Anti-Müllerian hormone (AMH) has become the ‘molecule of the moment’ in the field of reproductive endocrinology. Serum AMH concentrations show a progressive decline with female ageing. Indeed, it is valuable as a means of increasing understanding of ovarian pathophysiology and for guiding clinical management across a broad range of conditions.¹ AMH is an effective measure of ovarian reserve and it can predict ovarian response to controlled stimulation for women. The clinical pregnancy rate tends to increase as AMH increases.² The serum AMH level is an accurate marker of the ovarian early antral follicle number and offers a good diagnostic potency. In situations where accurate ultrasonographic data are not available, AMH could thus be used instead of the follicle count as a diagnostic criterion and incorporated as such in the Rotterdam definition of PCOS.³

[PRINCIPLE]

The AMH Rapid Test Cassette (Whole Blood/Serum/Plasma) is a semi-quantitative membrane based immunoassay for the detection of anti-Müllerian Hormone in whole blood, serum or plasma. In this test procedure, anti-AMH antibodies are immobilized in the test line region of the test. After specimen is added to the specimen well of the cassette, if AMH is present, it reacts with anti-AMH antibody coated colloid gold particles in the test. This mixture migrates chromatographically along the length of the test and interacts with the immobilized anti-AMH antibodies. If the specimen contains AMH, a colored line will appear in the test line region indicating a positive result. If the specimen does not contain AMH, a colored line will not appear in this region indicating a negative result. To serve as a procedural control, a colored line will always appear in the control line region, indicating that proper volume of specimen has been added and membrane wicking has occurred.

[REAGENTS]

The test contains anti-AMH antibody coated particles and anti-AMH antibodies coated on the membrane.

[PRECAUTIONS]

- For professional *in vitro* diagnostic use only. Do not use after the expiration date.
- Do not eat, drink or smoke in the area where the specimens or kits are handled.
- Do not use test if pouch is damaged
- Handle all specimens as if they contain infectious agents. Observe established precautions against microbiological hazards throughout testing and follow the standard procedures for proper disposal of specimens.
- Wear protective clothing such as laboratory coats, disposable gloves and eye protection when specimens are being tested.
- The used test should be discarded according to local regulations.
- Humidity and temperature can adversely affect results.

[STORAGE AND STABILITY]

Store as packaged in the sealed pouch at room temperature or refrigerated (2-30°C). The test is stable through the expiration date printed on the sealed pouch. The test must remain in the sealed pouch until use. **DO NOT FREEZE.** Do not use beyond the expiration date.

[SPECIMEN COLLECTION AND PREPARATION]

- The AMH Rapid Test Cassette can be performed using whole blood (from venipuncture or fingerstick), serum or plasma.
- To collect **Fingerstick Whole Blood** specimens: □ Wash the patient's hand with soap and warm water or clean with an alcohol swab. Allow to dry.
- Massage the hand without touching the puncture site by rubbing down the hand towards the fingertip of the middle or ring finger.
- Puncture the skin with a sterile lancet. Wipe away the first sign of blood.
- Gently rub the hand from wrist to palm to finger to form a rounded drop of blood over the puncture site.
- Add the Fingerstick Whole Blood specimen to the test by using **a capillary tube**.
- Touch the end of the capillary tube to the blood until filled to approximately 75 µL. Avoid air bubbles.
- Place the bulb onto the top end of the capillary tube, then squeeze the bulb to dispense the whole blood to the specimen well of the test cassette.
- Separate serum or plasma from blood as soon as possible to avoid hemolysis. Use only clear non-hemolyzed specimens.
- Testing should be performed immediately after the specimens have been collected. Do not leave the specimens at room temperature for prolonged periods. Serum and plasma specimens may be stored at 2-8°C for up to 3 days. For long term storage, specimens should be kept below -20°C. Whole blood collected by venipuncture should be stored at 2-8°C if the test is to be run within 1 day of collection. Do not freeze whole blood specimens. Whole blood collected by fingerstick should be tested immediately.

- Bring specimens to room temperature prior to testing. Frozen specimens must be completely thawed and mixed well prior to testing. Specimens should not be frozen and thawed repeatedly.
- If specimens are to be shipped, they should be packed in compliance with local regulations covering the transportation of etiologic agents.
- EDTA K2, Heparin sodium, Citrate sodium and Potassium Oxalate can be used as the anticoagulant for collecting the specimen.

[MATERIALS] Materials provided

- Test Cassettes
- Droppers
- Buffer
- Color Card
- Package insert

Materials required but not provided

- Specimen collection containers
- Centrifuge
- Timer

For fingerstick whole blood



- Heparinized capillary tubes and dispensing bulb

[DIRECTIONS FOR USE]

Allow the test, specimen, buffer and/or controls to reach room temperature (15-30°C) prior to testing.

- Bring the pouch to room temperature before opening it. Remove the test cassette from the sealed pouch and use it as soon as possible.
- Place the cassette on a clean and level surface. For **Serum or Plasma** specimen:
 - Hold the dropper vertically and transfer **2 drops of serum or plasma** (approximately 50 µL) to the specimen well (S), then add **1 drop of buffer** (approximately 40 µL), and start the timer. See illustration below.For **Venipuncture Whole Blood** specimen:
 - Hold the dropper vertically and transfer **3 drops of whole blood** (approximately 75 µL) to the specimen well (S), then add **1 drop of buffer** (approximately 40 µL), and start the timer. See illustration below.For **Fingerstick Whole Blood** specimen:
 - To use a capillary tube: Fill the capillary tube and **transfer approximately 75 µL of fingerstick whole blood** specimen to the specimen well (S) of test cassette, then add **1 drop of buffer** (approximately 40 µL) and start the timer. See illustration below.
- Wait for the colored line(s) to appear. **Read results at 10 minutes by comparing the T line intensity with provided color card.** Do not interpret the result after 20 minutes.

Note: It is suggested not to use the buffer, beyond 6 months after opening the vial.

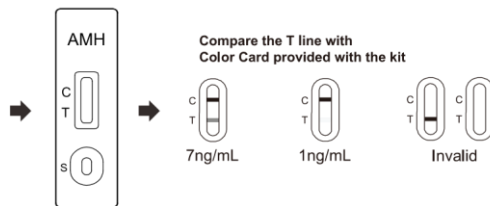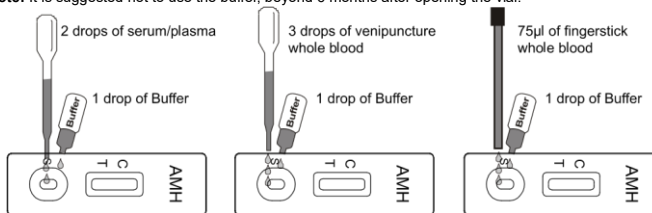

[INTERPRETATION OF RESULTS] (Please refer to the illustration and compare the T line intensity with “AMH Color Card” provided with the kit.)

AMH Level		Reference Range (ng/mL)
Normal		1-7
Abnormal	Indicative of lower ovarian reserve	<1
	Indicative of PCOS	>7

Normal: Two colored lines appear. One is in the control region (C) and another should be in the test region (T). **The line intensity in the test region (T) is darker than the 1ng/mL line and lighter than 7ng/mL line depicted on Color card provided with the kit.**

Abnormal :

Indicative of lower ovarian reserve : One colored line appears in the control region (C). No apparent colored line appears in the test region (T).

Indicative of PCOS: Two colored lines appear, one is in the control region (C) and another should be in the test region (T). **The line intensity in region (T) is darker than the 7ng/mL line depicted on the Color card.**

Note: Always compare the T line intensity with “AMH Color Card” and interpret results accordingly. **INVALID: Control line fails to appear.** Insufficient specimen volume or incorrect procedural techniques are the most likely reasons for control line failure. Review the procedure and repeat the test with a new test cassette. If the problem persists, discontinue using the test kit immediately and contact your local distributor.

[QUALITY CONTROL]

A procedural control is included in the test. A colored line appearing in the control region (C) is the internal procedural control. It confirms sufficient specimen volume and correct procedural technique. Control standards are not supplied with this kit; however, it is recommended that positive and negative controls be tested as a good laboratory practice to confirm the test procedure and to verify proper test performance.

[LIMITATIONS]

- The AMH Rapid Test Cassette (Whole Blood/Serum/Plasma) is for *in vitro* diagnostic use only. This test should be used for the detection of AMH in whole blood, serum or plasma specimen.
- The AMH Rapid Test Cassette (Whole Blood/Serum/Plasma) will only indicate the semi-quantitative level of AMH in the specimen and should not be used as the sole criteria for the diagnosis of ovarian functions.
- As with all diagnostic tests, all results must be interpreted together with other clinical information available to the physician.
- Physicians should diagnosis Polycystic ovary syndrome (PCOS) by combining a detailed history and the diagnostic criteria of the Rotterdam ESHRE/ASRM consensus for PCOS.

[EXPECTED VALUES]

The AMH Rapid Test Cassette (Whole Blood/Serum/Plasma) has been compared with CLIA demonstrating an overall accuracy of 97.1%.

[PERFORMANCE CHARACTERISTICS] Accuracy

The AMH Rapid Test Cassette (Whole Blood/Serum/Plasma) has been compared with CLIA. The following results are tabulated :

		Precision			Total Results
Method		CLIA			
AMH Rapid Test Cassette	Results	Lower (< 1ng/mL)	Normal (1-7ng/mL)	PCOS (>7ng/mL)	
	Lower (< 1ng/mL)	32	3	0	35
	Normal (1-7ng/mL)	2	149	1	152
	PCOS (>7ng/mL)	0	0	18	18
Total Results		34	152	19	205
Accuracy		94.1%	98.0%	94.7%	97.1%

Intra-Assay

Within-run precision has been determined by using 3 replicates of three specimens: 1ng/mL, 7ng/mL and 10ng/mL positive. The specimens were correctly identified >99% of the time. **Inter-Assay** Between-run precision has been determined by 3 independent assays on the same three specimens: 1ng/mL, 7ng/mL and 10ng/mL positive. Three different lots of the AMH Rapid Test Cassette (Whole Blood/Serum/Plasma) have been tested using 1ng/mL, 7ng/mL and 10ng/mL positive specimens. The specimens were correctly identified >99% of the time.

Cross-reactivity

The AMH Rapid Test Cassette (Whole Blood/Serum/Plasma) has been tested by inhibin-A, activin-A, LH, FSH, TGFβ-1 positive specimens. The results showed no cross-reactivity.

Substances

The AMH Rapid Test Cassette (Whole Blood/Serum/Plasma) has been tested for possible interference specimens. No interference was observed.

In addition, no interference was observed in specimens containing up to 300mg/dL Hemoglobin, 170 µmol/L ascorbic acid, 20 mg/dL acetaminophen, 65 mg/dL acetylsalicylic acid, 12 g/dL human serum Albumin and 1.4mmol/L Uric Acid.

[BIBLIOGRAPHY]

- A La Marca a,*, G Sighinolfi a, S Giulini a, M Traglia b, C Argento a, C Sala b, C Masciullo b, A Volpe a, D Toniolo b. Normal serum concentrations of anti-Müllerian hormone in women with regular menstrual cycles. a Mother–Infant Department, Institute of Obstetrics and Gynecology, University of Modena and Reggio Emilia, 41100 Modena, Italy; b Division of Genetics and Cell Biology, San Raffaele Scientific Institute, 4A2, Via Olgettina 58, 20132 Milano, Italy
- Sezai Sahmay & Mahmut Oncul & Abdullah Tuten & Abdullah Tok & Abdullah Serdar Acikgoz & Ismail Cepni. Anti-Müllerian hormone levels as a predictor of the pregnancy rate in women of advanced reproductive age. Received: 28 April 2014 / Accepted: 21 August 2014 / Published online: 4 September 2014# Springer Science+Business Media New York 2014
- P. Pigny, S. Jonard, Y. Robert, and D. Dewailly. Serum Anti-Müllerian Hormone as a Surrogate for Antral Follicle Count for Definition of the Polycystic Ovary Syndrome. Laboratoire de Biochimie et Hormonologie

(P.P.), Parc Eurasante ; Service de Gyne'cologie Endocrinienne et de Me'decine de la Reproduction (S.J., D.D.), Ho'pital Jeanne de Flandre; and Service de Radiologie, Ho'pital Jeanne de Flandre (Y.R.), Centre Hospitalier Re'gional Universitaire de Lille, F-59037 Lille, France

Index of Symbols

	Consult Instructions For Use		Tests per kit		Authorized Representative
	For in vitro diagnostic use only		Use by		Do not reuse
	Store between 2-30°C		Lot Number		Catalog #
	Do not use if package is damaged		Manufacturer		

 Hangzhou AllTest Biotech Co., Ltd.
#550, Yinhai Street
Hangzhou Economic & Technological Development Area
Hangzhou • 310018, P. R. China
www.alltests.com.cn




MedNet GmbH
Borkstrasse 10
48163 Muenster
Germany

Number: 146151500
Effective Date: 2019-11-29